




**PHYSICAL
ELECTRONICS**
A DIVISION OF ULVAC-PHI

SmartSoft-VersaProbe 3.5 for Windows 11

SmartSoft-VersaProbe 3.5 is the latest Windows 11 release of PHI data acquisition software with new features for the PHI VersaProbe Hybrid, PHI VersaProbe II, PHI VersaProbe III, and VersaProbe 4. SmartSoft-VersaProbe 3.5 is also compatible with Windows 10.

SmartSoft-VersaProbe Key Features

- New *SmartSoft* Navigator for novice users with preset acquisition parameters and enhanced flexibility for expert users
- New header bar for real-time status indicators in a single dashboard: loaded samples, data acquisition progress, vacuum pressure, ion guns, etc.
- Auto-tuning for monoatomic argon ion gun (Ar+) and gas cluster ion beam gun (GCIB)
- Intuitive queue editor for automation & throughput
- Software-guided intro photo alignment for precise point-and-click in both intro photo and scanning X-ray image (SXI)
- Simplified data acquisition setup for advanced techniques including ultraviolet photoelectron spectroscopy (UPS)
- Auto-tuning for ion gun, electron gun, and energy analyzer
- Multilevel system status logging and password protection for critical system parameters
- Greater control of automated chemical mapping and mosaic via SmartMosaic: sample height, rotation, and tilting
- Updated *SmartSoft-VersaProbe* Data Acquisition Software

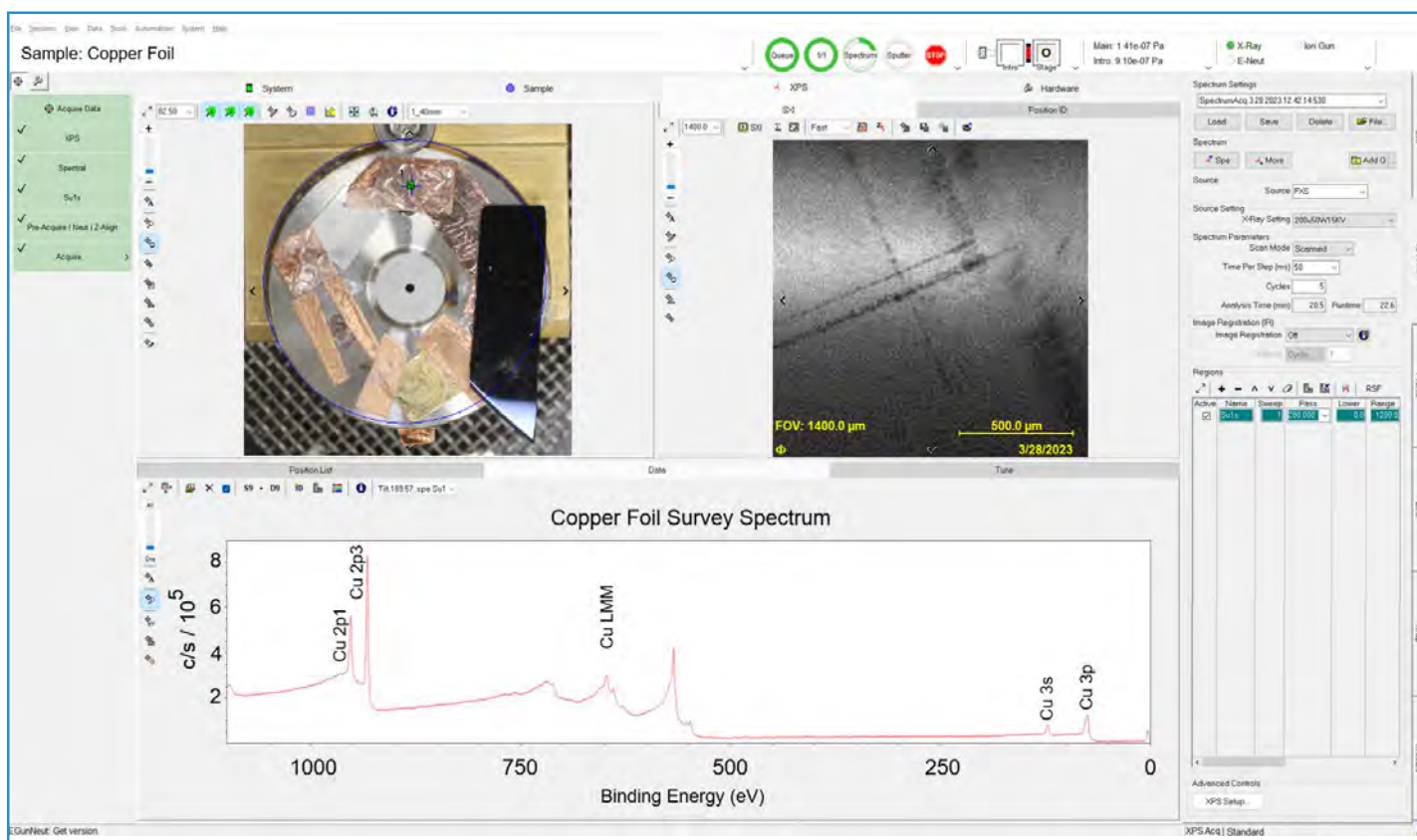

Updated SmartSoft-VersaProbe Software

**PHYSICAL
ELECTRONICS**
A DIVISION OF ULVAC-PHI

18725 Lake Drive East
Chanhassen, MN 55317, USA

www.phi.com
www.ulvac-phi.com

2500 Hagisono, Chigasaki
253-8522, Japan